

## On-Demand Cell Encapsulation Using On-Demand Droplet Generation and Impedance-based Detection

Tech ID: 21236 / UC Case 2011-162-0

### BRIEF DESCRIPTION

Current methods for cell encapsulation into droplets are inefficient and result in empty and unusable droplets. Researchers at UC Irvine have invented a novel microfluidic device capable of encapsulating cells at a very high efficiency. This tool can be used as a research tool for sorting and classifying cells for further downstream cell assays.

### FULL DESCRIPTION

Cell encapsulation into droplets is an integral part of droplet microfluidics. Improved cell encapsulation into the droplets is necessary as current methods for cell loading into droplets are inefficient, resulting in a majority of empty and unusable drops. Researchers at UC Irvine have invented a novel microfluidic system for generating droplets with highly efficient cell encapsulation.

The presence of objects such as cells in a fluid changes the electrical properties of the fluid. The invented system integrates impedance measurements with on-demand droplet generation process to enable selective generation of droplets when the presence of an object such as a cell is detected. This ensures that a high percentage of droplets that are generated actually contain the object and addresses the persistent problem of empty droplets. The impedance sensor not only determines object presence, but also can determine the size, shape of the object, and speed at which the object is moving.

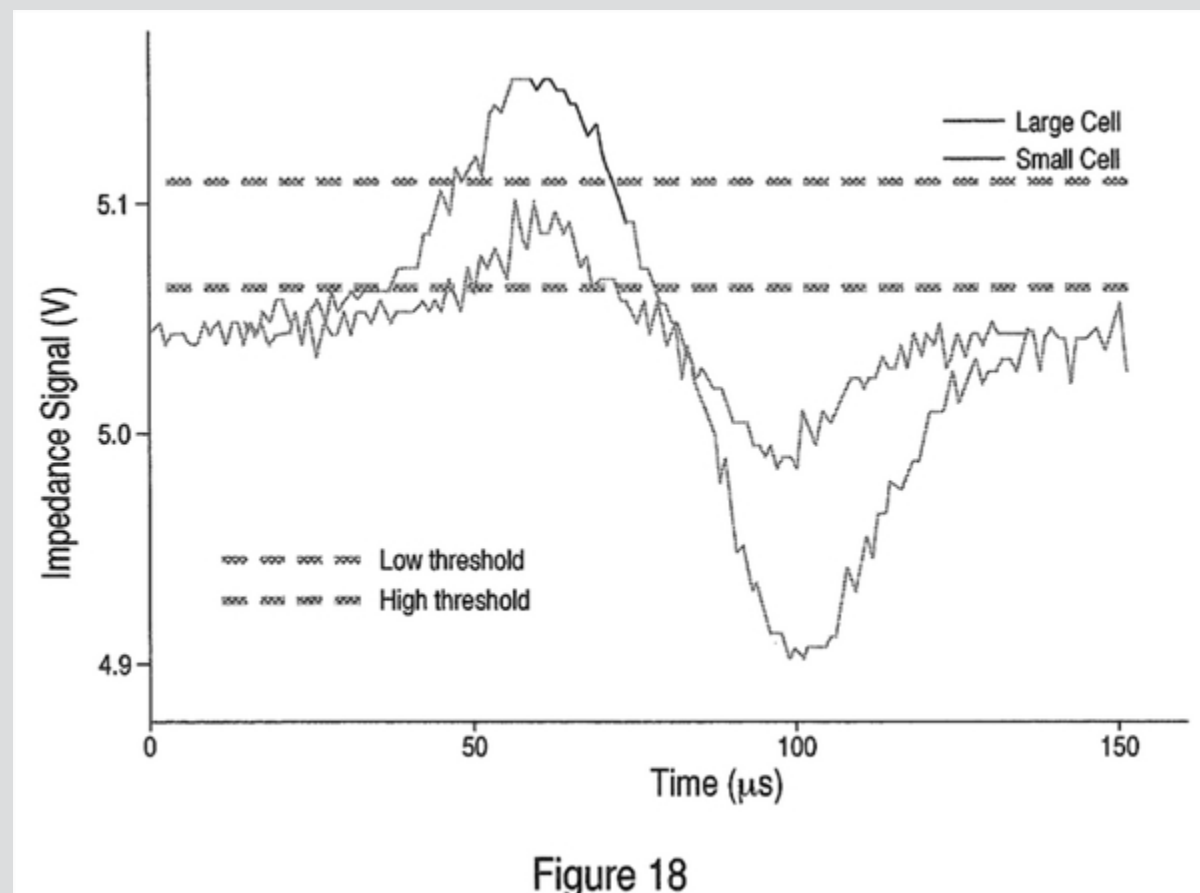


Figure 18 shows a graph of cell impedance data that is used for impedance signal processing. The present invention is able to distinguish object size and speed and use that information to generate the appropriate

### CONTACT

Alvin Viray  
aviray@uci.edu  
tel: 949-824-3104.

INTRODUCING  
**UC TechAlerts**  
New technology matches delivered to your email at your preferred schedule  
SEARCH ▶ SAVE SEARCH  
Learn More

### INVENTORS

- » Lee, Abraham P.
- » Lin, Robert
- » Lopez-Prieto, Javier

### OTHER INFORMATION

#### KEYWORDS

microfluidics, cell encapsulation, cell sorting, droplets

#### CATEGORIZED AS

- » **Biotechnology**
- » Other
- » **Medical**
- » Devices
- » Research Tools

sized droplet.

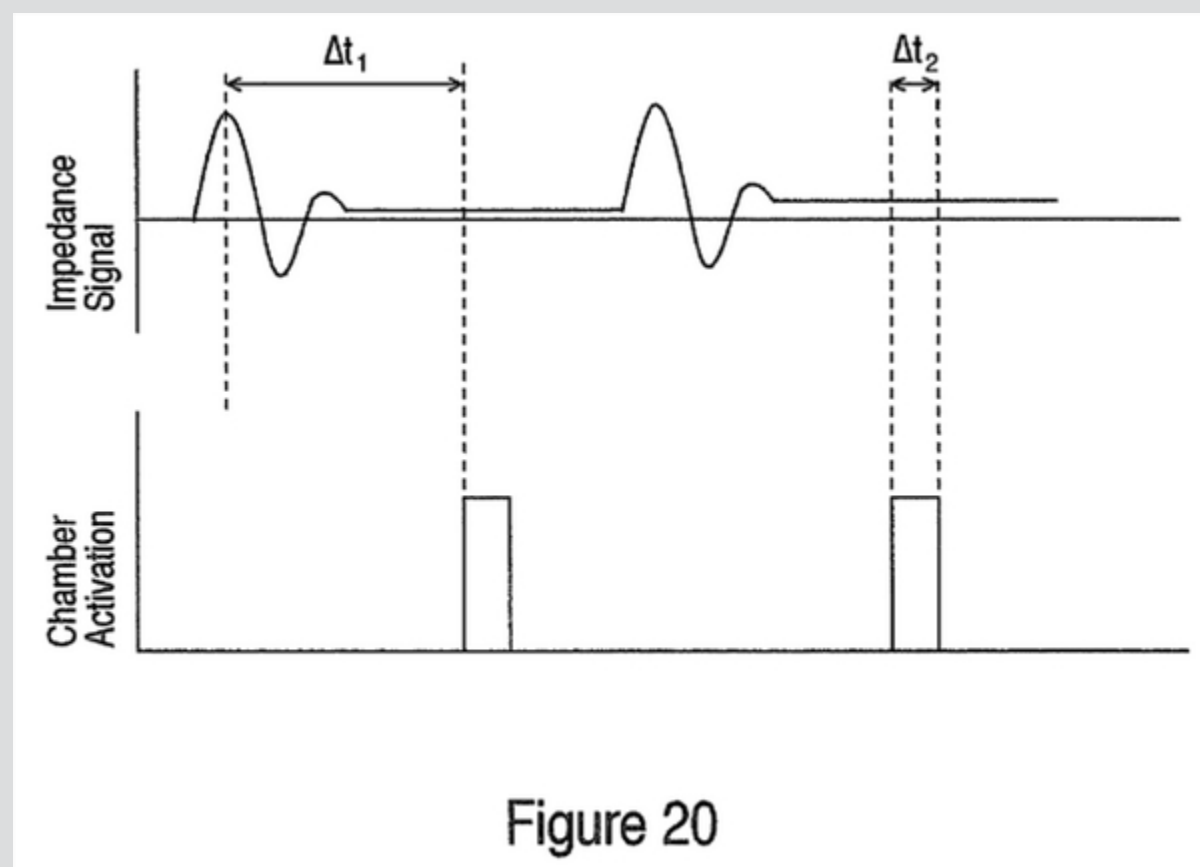


Figure 20

Figure 20 illustrates the encapsulation control scheme of the invention. Following impedance, or detection of cell presence, the chamber is activated for cell encapsulation.

The present invention has been shown to be able to accommodate densities of  $1 \times 10^6$  cells/mL.

## SUGGESTED USES

This microfluidic system can be used as a droplet encapsulation device. It can also be used for selective extraction of precise, digitized volumes out of a larger, continuous volume of fluid. Since the invention utilizes an impedance-based detector, this invention can also be used to select for cells with unique impedance signatures, for example, encapsulating white blood cells from whole blood.

## ADVANTAGES

The system senses changes in electrical properties of the fluid when objects such as cells are present, and generates droplets when cell presence is confirmed. This ensures a highly efficient process that can withstand a high density of cells. This process can distinguish between sizes and speed of object type and account for differences by changing the size of the produced droplet.

## PATENT STATUS

Country	Type	Number	Dated	Case
United States Of America	Issued Patent	9,999,886	06/19/2018	2011-162

## ADDITIONAL TECHNOLOGIES BY THESE INVENTORS

- ▶ New Microwell Plate Configurations to Increase Microwell Density
- ▶ Multi Layered Microfluidic Devices For In Vitro Large Scale Perfused Capillary Networks
- ▶ Controlled 'One-Cell-One-Bead' Encapsulation in Droplets
- ▶ Microfluidic Flow Transducer Based on the Measurement of Electrical Admittance
- ▶ Microfluidic device for multiplex diagnostics / Microfluidic devices and methods
- ▶ Microfluidic Device for Cell Separation Using Dielectrophoresis and/or Magnetohydrodynamics
- ▶ High throughput and precision cell sorting
- ▶ High-throughput Microfluidic Research Platform for Performing Versatile Single-Cell Molecular Timed-Release Assays within Droplets

» Screening

» **Nanotechnology**

» Other

» Tools and Devices

» **Research Tools**

» Antibodies

» Cell Lines

» Other

## RELATED CASES

2011-162-0

**UCI** Beall  
Applied Innovation

5270 California Avenue / Irvine, CA  
92697-7700 / Tel: 949.824.2683



© 2010 - 2018, The Regents of the University of  
California  
[Terms of use](#)  
[Privacy Notice](#)