UCI Beall Applied Innovation

Research Translation Group

Research Translation Group

Available Technologies

Contact Us

Request Information

Permalink

Controlled 'One-Cell-One-Bead' Encapsulation in Droplets

Tech ID: 28925 / UC Case 2017-814-0

BRIEF DESCRIPTION

The invention relates to a technique using a microfluidic device to perform co-encapsulation of samples in droplets and sorting of said droplets. This technique can achieve a one-one-one encapsulation efficiency of about 80% and can significantly improve the droplet sequencing and related applications in single cell genomics and proteomics.

.

FULL DESCRIPTION

Droplet-based single cell assays are based on the ability to encapsulate and confine single cells in individual droplets and enable high-efficiency genome wide expression profiling. Encapsulation of a single cell and bead reagent is of critical importance when conducting droplet-based analysis of single cells. Currently, single cell encapsulation in droplets is randomly performed and not well controlled. One-cell-one-bead encapsulation efficiency is limited by an uncontrolled random droplet encapsulation of cells and beads (0.1% or 1 in 1000 droplets). Such a low droplet encapsulation efficiency hinders the realization of high throughput droplet-based experiments. Hence there is a need for improved microfluidic device for encapsulation in single droplets and sorting thereof.

UCI researchers have developed a microfluidic droplet encapsulation method to efficiently encapsulate one-cell and one bead in a single droplet. By designing a microfluidic device with controllable fluid flow and a high force liquid-liquid interface, one-cell one-bead droplet encapsulation efficacy can be well controlled. The invention features a two-step technique to achieve high efficiency cell indexing coupled with down-stream sorting (Figure 1). In the first step, >30% cell encapsulation can be achieved by an interfacial shearing method utilizing laminar flows and high shear liquid-liquid interface at a microfluidic junction. A sorting module is incorporated as a second step to improve the encapsulation efficiency further, by up to 80% by removing empty droplets and/or droplets that do not contain the desired number of cells, beads or both.

Utilizing the invented technology, it has been shown that a significant increase in droplet encapsulation efficiency can be readily achieved. While existing techniques only provide random encapsulation governed by Poisson statistics, the improved encapsulation efficiency in this invention enhances the ability to perform single cell genomics and proteomics on a microfluidic platform.

Figure 1

CONTACT

Alvin Viray aviray@uci.edu tel: 949-824-3104.



INVENTORS

>> Lee, Abraham P.

OTHER INFORMATION

CATEGORIZED AS

» Medical

- » Diagnostics
- » New Chemical Entities, Drug Leads
- » Research Tools
- » Screening

» Research Tools

- Screening Assays
- » Sensors & Instrumentation
 - >> Process Control
 - » Scientific/Research

2017-814-0

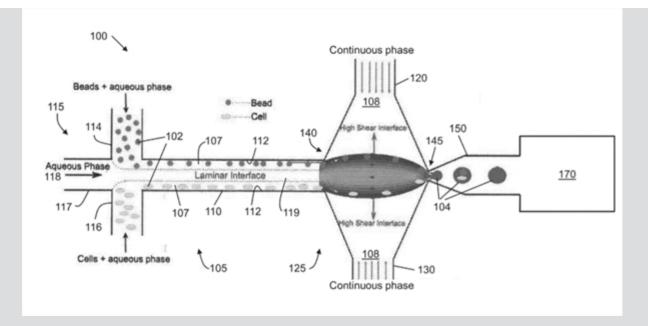


Figure 1 shows a schematic of high efficiency cell indexing in one-cell-one-bead droplets (1-1-1). Beads and cells introduced from upper and lower inlets self-assemble along the channel wall while moving toward high shear interfaces. At the droplet generation junction, both the beads and cells get pulled toward the high shear interface symmetrically from both the channel boundaries resulting in one-one-one encapsulation. A sorting module is incorporated downstream to sort out empty droplets from the bead-cell droplets.

STATE OF DEVELOPMENT

Inventors have completed the experimental and validation phase with a working prototype. As the next step, the inventors plan to optimize device parameters in hopes to commercialize or license the technology.

FEATURES/BENEFITS

Modular microfluidic technology that is capable of encapsulating one cell and one bead in one droplet

Increased droplet encapsulation efficiency

Improved droplet throughput for sequencing applications including single cell genomics and proteomics.

SUGGESTED USES

- >> Improving pharmacological assay throughput by efficiently generating droplets that encapsulate a single cell and bead
- >> Capture specific proteins and antigens from single cells using functionalized beads

PATENT STATUS

Country	Туре	Number	Dated	Case
United States Of America	Published Application	20200108393	04/09/2020	2017-814

Additional Patent Pending

ADDITIONAL TECHNOLOGIES BY THESE INVENTORS

- New Microwell Plate Configurations to Increase Microwell Density
- Multi Layered Microfluidic Devices For In Vitro Large Scale Perfused Capillary Networks
- Microfluidic Flow Transducer Based on the Measurement of Electrical Admittance
- Microfluidic device for multiplex diagnostics / Microfluidic devices and methods
- Microfluidic Device for Cell Separation Using Dielectrophoresis and/or Magnetohydrodynamics
- ▶ On-Demand Cell Encapsulation Using On-Demand Droplet Generation and Impedance-based Detection
- ► High throughput and precision cell sorting
- ► High-throughput Microfluidic Research Platform for Performing Versatile Single-Cell Molecular Timed-Release Assays within Droplets

UCI Beall Applied Innovation

5270 California Avenue / Irvine,CA 92697-7700 / Tel: 949.824.2683



© 2017 - 2020, The Regents of the University of California Terms of use Privacy Notice